

4 Ch Obese Rat Thorax Array *

MR Coils Made-to-Measure

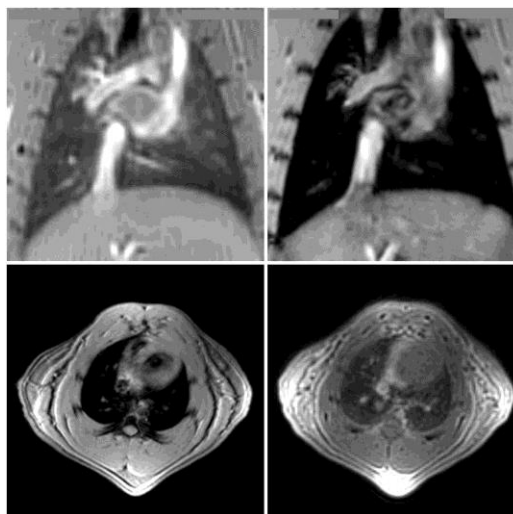


Image Courtesy Detlef Stiller, Boehringer Ingelheim Pharma GmbH & Co. KG Div. Research, Germany

- 4 channels for optimised image quality
- whole thorax coverage by volumetric array
- split coil housing supports
 - easy animal handling
 - reproducible animal positioning
 - fast animal change for maximum throughput
- integrated preamplifiers for optimum sensitivity
- large diameter for obese rats

Specifications

coil	rat thorax array coil
features	volumetric coil array 4 channels split housing for easy animal handling
B_0 field strength	1.5 T – 9.4 T
required equipment	actively decoupled volume transmit resonator with minimum inner diameter 72 mm