

PRODUCT DATASHEET

Non-Proton Torso Coils


Specifications

- volume transmit coil for human torso e.g. prostate, kidney, liver, heart (^{13}C) or lung (^{129}Xe)
- other nuclei on request
- 4 rung birdcage
- quadrature polarization
- transmit and transmit/receive
- open design to adapt bigger patients
- upper part detachable for patient access
- for 3 T and 7 T MR systems



transmit coil

add a dual tuned receive array for increased SNR

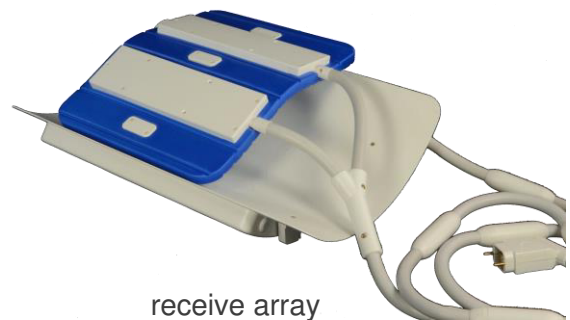

receive array

Regulatory requirements for medical devices will vary by country and MR system.

Please, contact us at info@rapidbiomed.de or info@rapidmri.com (USA, Canada, and South America) to determine approval status for products mentioned on this product sheet.